

### **OMNI**

# OMNI is with you Fight the virus

25g/m<sup>2</sup>

Thickened meltblown cloth

3layer

Triple filtering

GB-2010

International standard



Porous and breathable



Soft and comfortable



Low sensitivity and non-irritating



### Be careful!

## Inferior masks are more terrible than pollution!

Cheap masks use unknown raw materials and have no filtering and protective effects. The harmful substances in them, on the contrary, will harm the human body.



# Buying a mask is buying a meltblown cloth!

Masks that do not use meltblown cloth as filter material are all hooligans




25g/m<sup>2</sup>

Thickened meltblown cloth

0.5µm

Level protection

### 3 levels Protect for health!

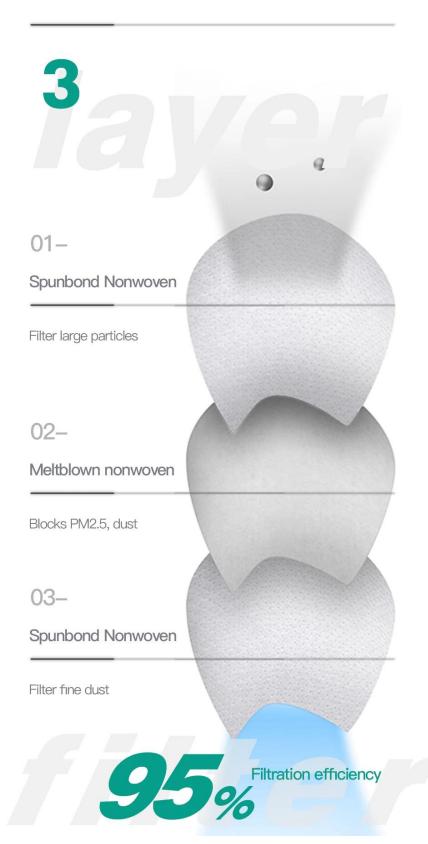

### SAFE

## Every breath you take needs to be supported safely



#### GB 18401-2010

Formaldehyde detection meets standards



#### GB 19083-2010

Microbiological and skin irritation tests meet standards

### Surrounded by softness

Non-woven inside and outside skin, as soft as undergarments




Nonwoven fabric

## Comfort 5le

Highly breathable and non-irritating

The surface is densely packed with air holes, comfortable and breathable, low sensitivity and no stimulation

## Three-dimensional protection

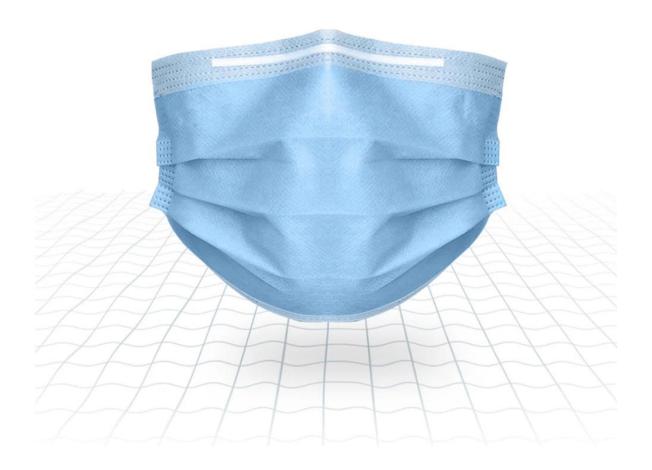

Length and width: 17.5x9.5cm

### Medium size and not tighten crowded

Suitable for most face shapes, comfortable fit, cotton elastic band, comfortable ears

## We all care about you for you

Detail

Plastic nose bridge

Fitted protection design

Wide cotton straps

### Information

Model: flat

Specifications: 17.5x 9.5cm

#### [Product performance structure and composition]

This product consists of a mask body, a nose clip, and a mask band. The cover body is composed of three layers of filter materials; the front and back layers are polypropylene spunbond non–woven fabrics, the middle layer is polypropylene melt–blown non–woven fabric, and the nose clip is metal core plastic. The mask band is made of elastic polyester spandex. For non–sterile products,—for single use.

#### [Product scope and use]

It is suitable for covering the user's mouth, nose and jaw. It is used in ordinary medical environment to wear and block the mouth and nasal cavity from exhaling or spraying pollutants.

#### [Instructions for use]

Flatten the mask with the dark side facing outwards and the nose clip facing up. Wear it on the nose and mouth and hang the earband on the ear.

### [Precautions, warnings, and informative instructions]

- Please use within the validity period.
- 2. Only for one-time use, please destroy and dispose according to regulations after use.
- 3. The packaging is damaged and it is forbidden to use it.
- 4 Use with caution to those who are allergic to nonwovens.
- 5 Please read the instructions before use.

### [Special storage conditions and methods]

Store in a relative humidity not exceeding 80%, non-corrosive gas and well-ventilated room, avoid high temperature.

